




309 Laurelwood Rd. Suite 20  
Santa Clara, CA 95054  
sales@sstgroup-inc.com

Sales Support: 800.944.6281  
Outside US: 408.350.3450  
Fax: 408.350.3100

## DMC - EZ 575

The DMC-EZ 575 DICOM Publishing system is perfect for clinics, imaging centers, and hospitals. It allows users to create DICOM CDs and DVDs for distribution to patients and referring physicians. The DMC-EZ systems satisfy the requirements of all radiology and PACS systems while giving users an intuitive system that is robust and easy to use.

## DICOM Image Distribution

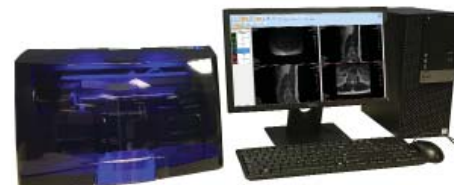

Speed Capability: 30 Full CDs per hour

### DMC-EZ Medical Software

- Intuitive browser-based user interface.
- Compatible with large format studies (TOMO, Cardio and Ultrasound)
- Creates IHE DICOM-compliant CDs and DVDs.
- On-demand CD/DVD creation for DICOM studies (Query/Retrieve).
- On-demand cloud delivery of DICOM studies (HITECH compliant) - optional module.
- Automatic disc creation using powerful rule-based scheduling utilizing multiple AE Titles.
- Includes multiple viewers for disc distribution. Windows viewer supports XP through Windows 10.
- Includes diagnostic reports on CD via HL7 feed, SR, Mitra Broker. We also integrate with major RIS/PACS OEM to obtain diagnostics reports. Contact sales for the options.
- View complete report on disc with images and hyperlinks.
- Flexible disc labeling options.
- Study disc anonymization.
- View thumbnail images prior to disc creation.
- Built-in web server and web viewer for alternative distribution of images and reports.
- Optional REST API allows burner to be integrated into your workflow.
- Image distribution via USB flash for easy use by referring physicians.
- Automated or Manual study disc encryption (HIPAA/HITECH 256 AES/FIPS)
- Manual and Automated Rules-based DICOM routing.
- Robot supplies status monitor
- All activity recorded in audit logs (HIPAA Compliance)
- Back-up configuration and HIPAA logs to USB, CD/DVD or network location to allow easy restore.
- Import data from CDs and DVDs your facility receives (EzDicom Sender).
- Multiple language support - interface and viewer (English, French, Italian, Portuguese, Spanish).
- User login supports active directory.
- **Patented technology.**

### Primera Bravo 4202 Disc Publishing System

High performance dual CD/DVD recorder inkjet desktop disc publishing system with high speed inkjet printer, automated robotic print and CD/DVD recording. Automatically burns and prints up to 100 discs unattended. (Depending on configuration)

### DMC-EZ DICOM Server and Software

High performance CPU contains the burning and printing software, DMC-EZ application and image cache. DMC-EZ PC running Windows® 7 Pro or Windows® 10 Pro controls Primera 4202 and connects to DICOM network or DICOM workstation. Includes display, keyboard and mouse.

### Color Inkjet Printing on CD/DVD Media

Discs surface printed at 4800 DPI. Inkjet media printing is now enhanced by the use of SST Group's MediScribe InstaDry™, Water Repel media which helps prevent degradation of the printed media surface when exposed to water or liquids.

### Fully Automated Disc Creation

Includes two CD/DVD recorders, two 50 disc input and one 100 disc output bin.

### Primera Bravo 4202

Power: 100VAC - 240VAC ±10%, 50/60Hz  
Operating Temp: Temperature: 60° - 86°F (16°C - 30°C)  
Physical Specifications: 17.8" W x 11" H x 17.5" D (452mm W x 279mm H x 445mm D)  
26 Lbs. (11.8 Kg)

### DMC-EZ PC

Power: 100 - 240VAC ±10%, 50/60Hz  
Operating Temp: Temperature: 60° - 86°F (16° - 30°C)  
Network: 100/1000 Ethernet interface  
Physical Specifications:  
6.1" W x 13.8" H x 10.8" D (154mm W x 350mm H x 274mm D)  
17.62Lbs. (8 Kg)