



Your DICOM distribution experts.

309 Laurelwood Rd. Suite 20
Santa Clara, CA 95054
sales@sstgroup-inc.com


Sales Support: 800.944.6281
Outside US: 408.350.3450
Fax: 408.350.3100

www.sstgroup-inc.com

DMC - EZ 575

DICOM Image Distribution

The DMC-EZ 575 DICOM Publishing system is perfect for clinics, imaging centers, and hospitals. It allows users to create DICOM CDs / DVDs / USB and Virtual Media (Cloud) for distribution to patients and referring physicians. The DMC-EZ systems satisfy the requirements of all radiology and PACs systems while giving users an intuitive system that is robust and easy to use.

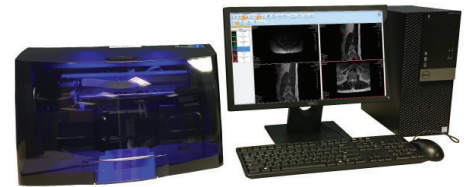

Speed Capability: 30 Full CDs per hour

DMC-EZ Medical Software

- Intuitive browser-based user interface.
- Compatible with large format studies (TOMO, Cardio and Ultrasound)
- Creates IHE DICOM compliant CDs and DVDs.
- On-demand CD/DVD creation for DICOM studies (Query/Retrieve).
- **On-demand cloud delivery of DICOM studies (HITECH compliant) - optional software. (Zero Trust Image Sharing)**
- **New - Create IHE DICOM compliant USB with labeling - optional hardware and software.**
- Automatic disc creation using powerful rule-based scheduling utilizing multiple AE Titles.
- Includes multiple viewers for disc distribution. Windows viewer supports XP through Windows 11.
- Includes diagnostic reports on CD via HL7 feed, SR, Mitra Broker. We also integrate with major RIS/PACS OEM to obtain diagnostics reports. Contact sales for the options.
- View complete report on disc with images and hyper-links.
- Flexible disc labeling options.
- Study disc anonymization.
- View thumbnail images prior to disc creation.
- Built-in web server and web viewer for alternative distribution of images and reports.
- Optional REST API allows burner to be integrated into your workflow.
- Automated or Manual study disc encryption (HIPAA/HITECH 256 AES/FIPS)
- Manual and Automated Rules-based DICOM routing.
- Robot supplies status monitor
- All activity recorded in audit logs (HIPAA Compliance)
- Back-up configuration and HIPAA logs to USB, CD/DVD or network location to allow easy restore.
- Import data from CDs, DVDs and USB your facility receives (EzDicom Sender).
- Multiple language support - interface and viewer (English, French, Italian, Portuguese, Spanish).
- User login supports active directory.
- **Patented technology.**

Primera Bravo 4202 Disc Publishing System

High performance dual CD/DVD recorder inkjet desktop disc publishing system with high speed inkjet printer, automated robotic print and CD/DVD recording. Automatically burns and prints up to 100 discs unattended. (Depending on configuration)

DMC-EZ DICOM Server and Software

High performance CPU contains the burning and printing software, DMC-EZ application and image cache. DMC-EZ PC running Windows® 10 Pro or 11 Pro controls Primera 4202 and connects to DICOM network or DICOM workstation. Includes display, keyboard and mouse.

Color Inkjet Printing on CD/DVD Media

Discs surface printed at 4800 DPI. Inkjet media printing is now enhanced by the use of SST Group's MediScribe InstaDry™, Water Repel media which helps prevent degradation of the printed media surface when exposed to water or liquids.

Fully Automated Disc Creation

Includes two CD/DVD recorders, two 50 disc input and one 100 disc output bin.

Primera Bravo 4202

Power: 100VAC - 240VAC $\pm 10\%$, 50/60Hz
Operating Temp: Temperature: 60° - 86°F
Physical Specifications: 17.8" W x 11" H x 17.5" D, 26 Lbs.

DMC-EZ PC

Power: 100 - 240VAC $\pm 10\%$, 50/60Hz
Operating Temp: Temperature: 60° - 86°F
Network: 100/1000 Ethernet interface
Physical Specifications: 6.1" W x 12.8" H x 11.5" D 14.59 lbs.