



Your DICOM distribution experts.

309 Laurelwood Rd. Suite 20
Santa Clara, CA 95054
sales@sstgroup-inc.com

Sales Support: 800.944.6281
Outside US: 408.350.3450
Fax: 408.350.3100

www.sstgroup-inc.com

DICOM Image Distribution

DMC - EZ DICOM Distribution Station


CD/DVD, USB & Cloud Delivery of Patient Images

EZ to implement and EZ to use, the DMC-EZ DICOM Distribution Station is the ultimate control center for supplying medical images to patients and physicians.

One simple interface allows you to select your patient images and send to stake holders via the method they prefer: optical media (CD/DVD), USB, cloud (via email).

A few clicks and your images are ready for patients and physicians whether they are local or off-site.

Simple. Powerful. Secure.

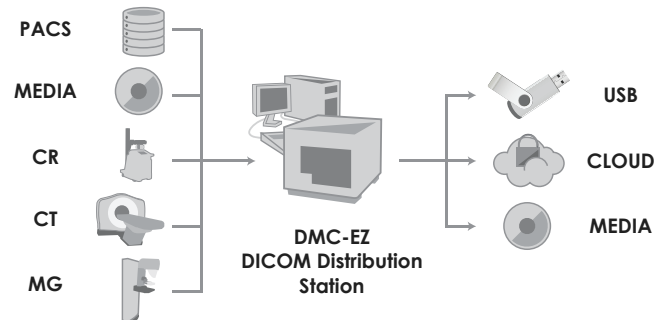

DMC-EZ DICOM Medical Software

- Intuitive browser-based user interface
- Compatible with large format studies (TOMO, Cardio and Ultrasound)
- Create IHE DICOM-compliant CD/DVD/USB with labeling
- On-demand CD/DVD/USB creation for DICOM studies (Query/Retrieve)
- On-demand cloud delivery of DICOM studies (HITECH compliant)
- Automatic CD/DVD/USB creation using powerful rule-based scheduling utilizing multiple AE Titles
- Includes multiple viewers for CD/DVD/USB distribution, Windows® viewer supports XP through Windows 11
- Includes diagnostic reports on CD/DVD/USB/virtual (cloud) via HL7 feed, SR, Mitra Broker and integrates with major RIS/PACS OEM to obtain diagnostics reports. (Contact sales for options.)
- View complete report on disc with images and hyperlinks
- Flexible disc labeling options
- Study CD/DVD/virtual/USB anonymization
- Built-in web server and web viewer for alternative distribution of images and reports
- Automated or Manual study encryption (HIPAA/HITECH 256 AES/FIPS)
- Robot supply status monitor
- All activity recorded in audit logs (HIPAA Compliance)
- Back-up configuration and HIPAA logs to USB, CD/DVD or network location for easy restore
- Import DICOM data from CD, DVD, USB or HDD (EzDicom Sender)
- Multiple language support for interface and viewer (English, French, Italian, Portuguese, Spanish)
- User login supports active directory
- **Patented technology**



DMC-EZ Server and Software

High performance CPU running Windows® 10 Pro or 11 Pro contains the burning and printing software which controls Epson Discproducer, DMC-EZ application and image cache. Connects to DICOM network and/or DICOM workstation(s).

Epson PP100 iii Discproducer

High performance dual CD/DVD recorder inkjet desktop disc publishing system with vibrant six-color inkjet printer and precision AcuGrip™ technology, automated robotic print and CD/DVD recording.

USB Components

4-port USB hub, label printer and 50-8GB SST branded USB drives and b/w printer ribbon.

Virtual Media Distribution

Simple, fast and secure method to get DICOM images to your patients, referring physicians and other facilities.

There are **no per study fees, no special sign-up or accounts to setup.**

(Virtual media supports Windows®7, 8, 10, 11)